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
Thank you for purchasing your beyondRED® PRO.

OVERVIEW / SCOPE OF DELIVERY

Make sure that all elements shown in the scope of delivery are included.

(1 – 10) Have fun with your device!

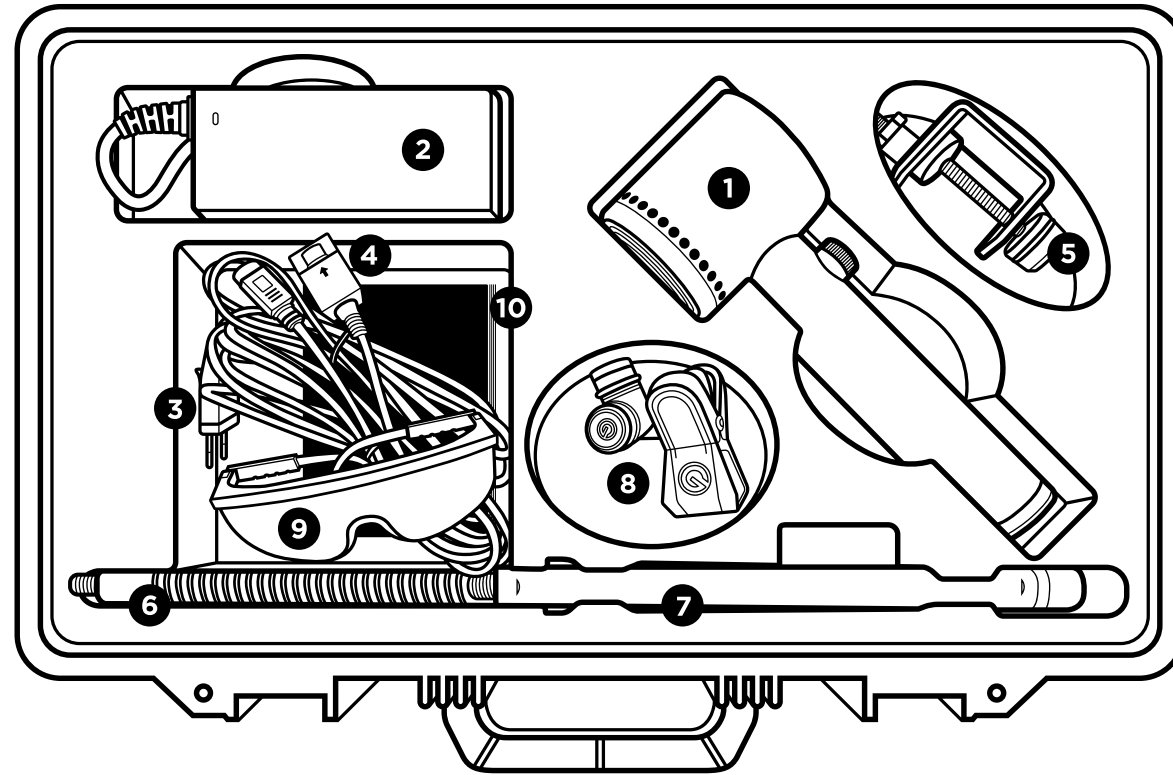
- 1 beyondRED® PRO
- 2 Power supply
- 3 Mains cable
- 4 Cable extension
- 5 Table clamp
- 6 Gooseneck
- 7 Extension rod
- 8 Device clamp
- 9 Safety goggles
- 10 Application recommendations


YOUR FIRST STEPS WITH THE BEYONDRED® PRO

1. POWER

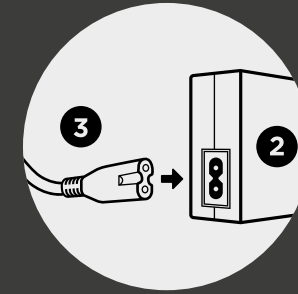
Connect the power supply

a) Remove the power supply 2 and connect it to the power cable 3.

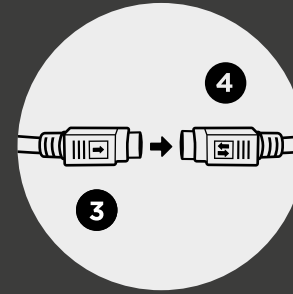b) As an extension, connect the mains cable 3 to the cable extension 4.
c) When connecting, pay attention to the directional arrows.


2. CHARGING

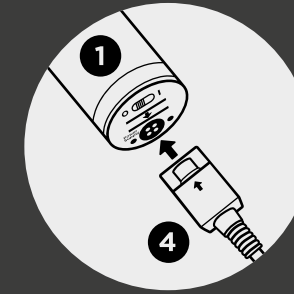
Charge the device fully before first use

a) Connect the cable extension 4 to the device 1. Pay attention to the correct position of the arrows.
b) The device should be charged for at least 2 hours.


3. BRACKET

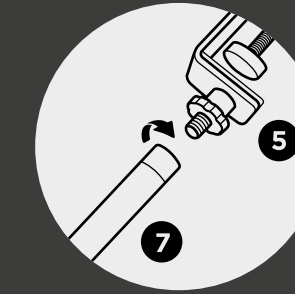a) Screw the table clamp 5 to the extension rod 7.

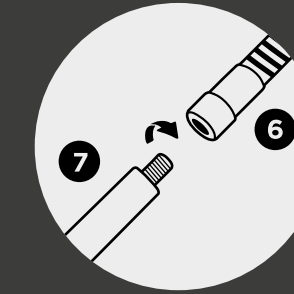
b) Screw the extension rod 7 to the gooseneck 6.

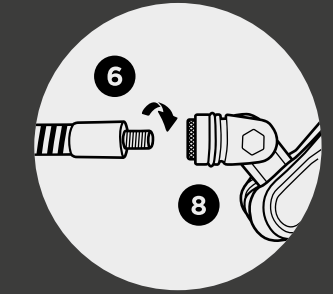
c) Connect the gooseneck 6 to the device clamp 8.

4. USE

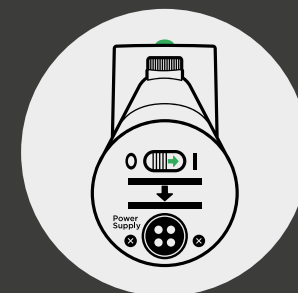
Now you are ready for treatment

a) Switch on the beyondRED® PRO using the power switch on the underside.
b) The status LED lights up green continuously.

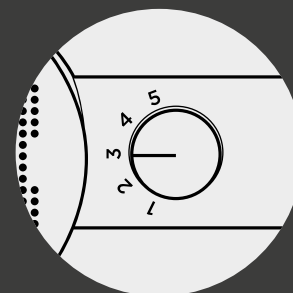
c) Turn the knob to set the desired intensity. Our application recommendations will help you with this 10.

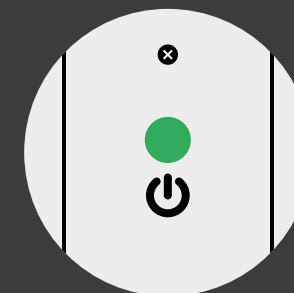d) Press the green treatment start button to start the treatment. The status LED flashes green during the entire treatment time.
e) The integrated fan starts up during treatment. This is used to cool the high-power diodes.
f) The beyondRED® PRO is equipped with a timer. The treatment ends automatically after approx. 5 minutes.

5. END

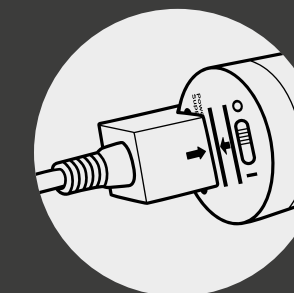End of treatment

a) After the treatment is over, you can switch off the device. To do this, set the power switch on the underside of the device to 0.

Other languages